

EWSR1 Break Apart FISH Probe Kit

Introduction

The EWSR1 Break Apart FISH Probe Kit is designed to detect rearrangements in the human *EWSR1* gene located on chromosome band 22q12.2. In addition to revealing breaks, which can lead to translocation of parts of the gene, inversion, or its fusion to other genes, the probe set can also be used to identify other *EWSR1* aberrations such as deletions or amplifications. Rearrangements and abnormal expression of the *EWSR1* gene – also known as *EWS* or *bK984G1.4* – have been observed in Ewing’s sarcoma and in numerous neuroectodermal and other tumor types.


Intended Use

To detect rearrangements in the human *EWSR1* gene located on chromosome band 22q12.2.

Cont.	Color
-------	-------

EWSR1 5' FISH Probe	CytoOrange CytoGreen
EWSR1 3' FISH Probe	

Probe Design

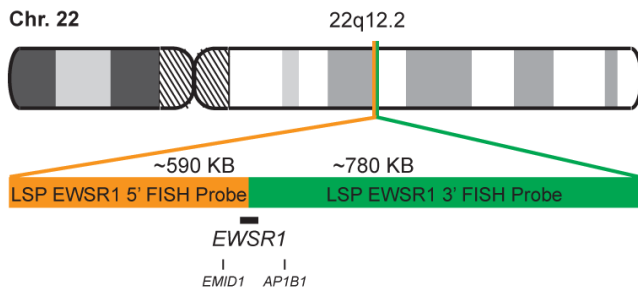

LSP EWSR1 5' FISH Probe covers the 5' (start) portion of the *EWSR1* gene and some adjacent genomic sequences. LSP EWSR1 3' FISH Probe covers the 3' (end) part as well as sequences downstream of the gene. The two probes are flanking sequences across the *EWSR1* gene in which variable breakpoints have been observed.

Not to Scale

Cat. No.	Volume
----------	--------

CT-PAC054-10-OG	10 Tests (100 µL)
-----------------	-------------------


Signal Pattern Interpretation

Normal Pattern	Abnormal Pattern
----------------	------------------

2OG*	Other Patterns
------	----------------

*Overlapping orange and green signals can appear as yellow.

- Ludwig JA. *Curr Opin Oncol.* 20(4):412-8 (2008).
- Erkizan HV, et al. *Clin Cancer Res.* 16(16):4077-83 (2010).
- Romeo S & Dei Tos AP. *Virchows Arch.* 456(2):219-34 (2010).
- Sohn EJ, et al. *Cancer Res.* 70(3):1154-63 (2010).
- Tanas MR, et al. *Mod Pathol.* 23(1):93-7 (2010).

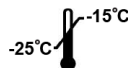

CytoTest Inc.
9430 Key West Ave., Suite 210
Rockville, MD 20850, USA

* CE IVD only available in certain countries. All other countries are either ASR or RUO. Please contact your local dealer or our headquarters for more information.